


Adapter Rings, Coverslips and Culture Dishes

Accessories and Replacement Parts for Warner Micro-incubation Chambers

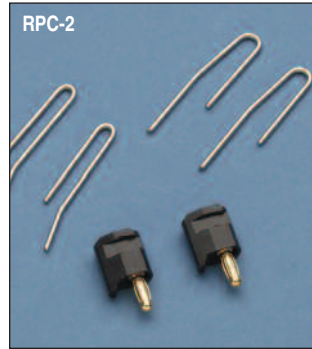

Order #	Model	Product
W3 64-0350	AR-1	Adapter Ring for DH-35 and Willco Dish
W3 64-1455	AR-2	Adapter Ring for AttoFluor® cell chamber in DH-35i/QE-1
W3 64-0351	AR-4	Adapter Ring Kit, set of 6, for Willco, Corning, Nunc, WPI, MatTek and Falcon Dishes
W3 64-1507	RPC-2	Replacement Clamps for DH-35i/DH-40i/QE-1
W3 64-1500	CS-40R	Cover Glass, 40 mm round, #1.0 thickness, 60 pc.

Heated Top Coverslip for DH-35i and DH-40i Chambers



Developed for use in our DH-35i and DH-40i Micro-incubation chambers, this heated coverslip assembly replaces the 40 mm coverslip normally supplied with these chambers. Once installed, the HCS-1 eliminates condensation which occurs when these chambers are used in the closed mode.

The heated coverslip derives its power from the heated base of the DH-35i or DH-40i chamber and does not require an additional temperature controller.

Order #	Model	Product
W3 64-1543	HCS-1	Heated Top Coverslip
W3 64-1544	HCS-RPL	Replacement ITO Coverslip

Glass Bottom Cell Culture Dishes (35 and 50 mm), #1.5 Thickness (0.17 mm)



- Highest quality construction, sterilized
- #1.5 Cover glass for microscope imaging
- Glass fully covers bottom
- Compatible with all Warner Dish Heaters

Glass Bottom Cell Culture Dishes

High quality polystyrene dishes with No. 1.5 cover glass bottoms are available in 35 mm and 50 mm diameters, available both in pouch packs (20/pouch) and individual blister packs (120/case). All dishes are sterilized before packaging.

50 mm dishes are 7 mm in height.

Order #	Model	Product
Pouch-Pack Dishes, 20 sterilized dishes/pouch		
W3 64-0757	D3512P	35 x 9.9 mm, 12 mm Well Diameter
W3 64-0758	D3522P	35 x 9.9 mm, 22 mm Well Diameter
W3 64-0759	D5030P	50 x 7.2 mm, 30 mm Well Diameter
W3 64-0760	D5040P	50 x 7.2 mm, 40 mm Well Diameter
Blister-Pack Dishes, individual packaged sterile dishes, case of 120		
W3 64-0761	D3512B	35 x 9.9 mm, 12 mm Well Diameter
W3 64-0762	D3522B	35 x 9.9 mm, 22 mm Well Diameter
W3 64-0763	D5030B	50 x 7.2 mm, 30 mm Well Diameter
W3 64-0764	D5040B	50 x 7.2 mm, 40 mm Well Diameter

Note: 35 mm dishes supplied with triple vented lids for optimum ventilation. 50 mm dishes have low side walls to facilitate microinjection.