

chamber accessories

Electrode Kits and Magnetic Clamps

Chambers

chamber accessories

Replacement parts and accessories for Series 20 chambers

Agar Bridge Reference Electrode Kit

Developed for use in low volume recording chambers. The ABR-1 Agar Bridge Reference Electrode Kit consists of a 1 mm diameter silver wire electrode mounted in a modified pipette tip that can be filled with a 3 M KCl agar mixture. An included magnetic clamp with flexible holder for the pipette tip is perfect for positioning the electrode at any angle in the chamber. In addition, the pipette tips can be cut to length or shaped using a heat gun. The silver wire electrode comes with 1 meter of flexible wire terminated in a 2 mm pin. Includes a 2 mm to 1 mm pin adapter for use with 1 mm jacks.



ABR-1 Kit Components




ABR-1 Shown mounted to a QE-1 Heated Base

Order #	Model	Product
W4 64-1426	ABR-1	Agar Bridge Reference Electrode Kit
W4 64-1428	LPE-10	Replacement Pipette Tips, pkg. of 10
W4 64-1323	PJ2-5	2 mm jack (5 pack)
W4 64-0360	MAG-3	Magnetic Clamp

Mini Magnetic Clamps

These miniature magnetic clamps offer a quick and easy method of positioning perfusion and suction tubes in chambers equipped with a magnetic platform. Sized to fit 20, 18, and 16 gauge tubes, they can also be used to hold PE tubing for applications requiring metal free solution delivery. Height: 16.6 mm, Ø: 6.4 mm.

Order #	Model	Product
W4 64-1552	MAG-5	Mini Mag for 20 gauge tubing
W4 64-1553	MAG-6	Mini Mag for 18 gauge tubing
W4 64-1554	MAG-7	Mini Mag for 16 gauge tubing

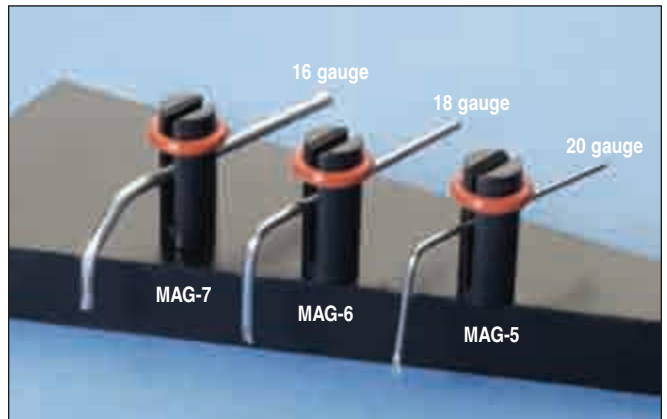

MAG-7 is for Suction, MAG-6 is for Perfusion, MAG-5 is for Skinny Suction Tubes

Suction/Perfusion Tubes

	A	
	A	B
ST-3	28.5 mm	6.5 mm
ST-3R	28.5 mm	6.5 mm
ST-QE1	40.0 mm ±	13.0 mm ±
PT-QE1	39.5 mm ±	12.7 mm ±