

Suction Tube Upgrade Kit

The STU-1 kit includes a mechanical clamp with two slits that hold on to small diameter tubing such as a suction tubes.

As shown in figure 1 each kit comes with a set of magnetic stainless steel mounting plates. The set can be attached to the platform clamps by using the double sided adhesive tabs.

Figure 2 shows how to align the plate on platform clamp, be sure to clean the clamp before mounting magnetic plate.

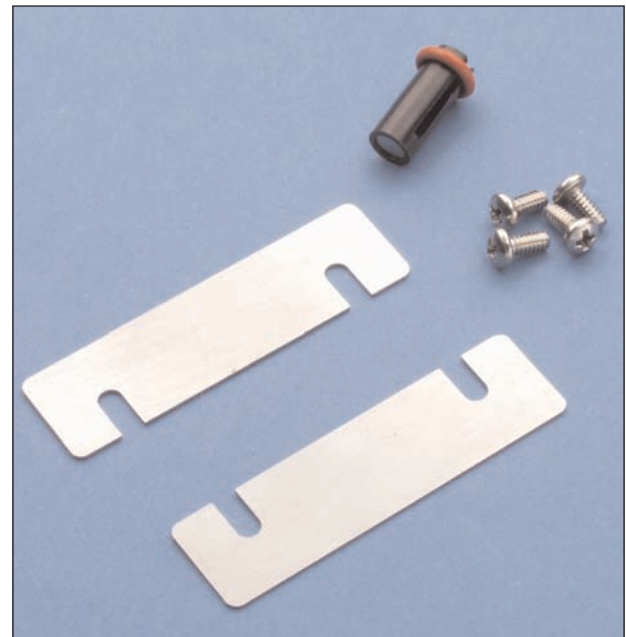


Figure 1. Magnetic Plates

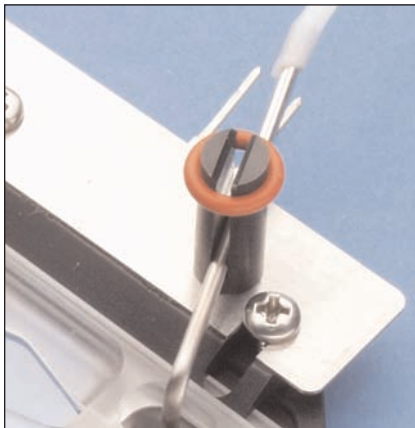


Figure 3. MAG-7 Suction Tube Holder

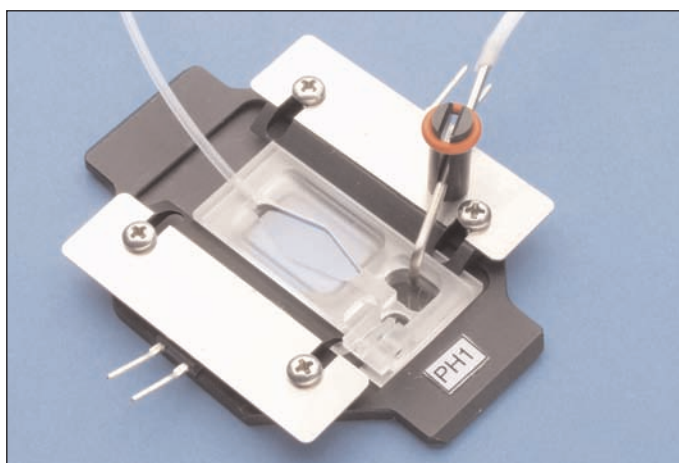


Figure 4. STU-1 Mounted to a P-1 Platform



Figure 2. Mounting Magnetic Plate using double sided adhesive tape